

# **China electric wheelchairs** **manufacturer**








# TRAVEL BUTLER SMART ELECTRIC WHEELCHAIR



## PRODUCT ADVANTAGES

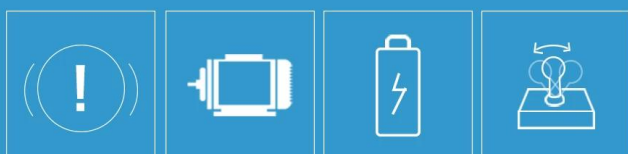

EABS Brake system

500W dual motor

High battery life

360° joystick



Hand/electric mode

Easy to fold

Aluminum alloy frame

Anti-backover design

## PRODUCT PARAMETER

Parameter	Category
Series/Model No.	Alphard Series/74501
Frame material	Aluminum alloy
Unfolded Dimension	1000*610*930mm
Folded Dimension	660*380*770mmmm
Motor	DC 24V、250W
Battery	Lithium battery DC 24V 13 AH
Weight	24KG
Loading Capacity	120KG
Front wheel specifications	8"*2" (200*50mm) PU tire
Rear wheel specifications	12"*2" (305*60mm) pneumatic tire
Max Speed	6KM/H
Driving Range	15KM

## PRODUCT STRUCTURE



Motor 1

Motor 2

500W dual motor,  
DC motor /strong climbing ability



LG cell lithium batteries  
travel distance about 15km

## PRODUCT DETAILS



Indicator light

Power switch

Horn

Power display

Speed adjustment

**360° CONTROLLER**

Wheelchair stop immediately when the hand leaves the controller.

### FOLDING DESIGN

Lightweight folding portable design  
Can be placed in the trunk, Similar to suitcase size.



### LEFT OR RIGHT HAND ADJUSTABLE

Can be adjusted as needed, the default is installed on the right



### THE ARMREST CAN BE RAISED

Handrails can swing up and down  
Easy to hold the patient up and down wheelchair



### HAND ELECTRIC SWITCH

Easily switch between manual and electric, multiple ways to travel

